


la belle CL 210



la belle CL 210




- BTE-module for eyeglasses
- module cold bendable
- digital hearing system
- natural sound transfer via ex-speaker system *
- 3 programmes selected by push button
- 2 channel WDRC
- Dynamic Contrast Detection
- 12-band gain adjustment
- MPO
- adaptive feedback cancellation (until -17dB)
- programme tones and battery tones
- ideally suited to open ear application
- takes 312 type batteries
- replaceable cerumen defenser filter
- tube warm bendable
- colours: gold-brown, silver-black, anthacite-black
- mono, stereo, CROS, BiCROS
- *ex-speaker = speaker in the ear canal

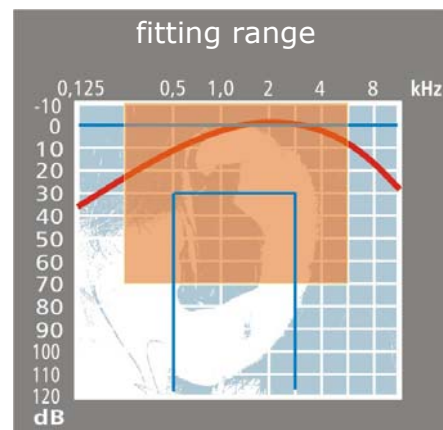

Technical data

| | | DIN EN 60118-0 | DIN EN 60118-7 |
|------------------------|------|----------------|----------------|
| OSPL90 max. | (db) | 119 | 110 |
| Gain 50dB(FOG) max. | (dB) | 49 | 42 |
| Frequency range | (Hz) | 200-6000 | 200-6500 |
| Equivalent input noise | (dB) | 22 | 22 |
| Operating current | (mA) | 1,2 | 1,2 |

please turn over

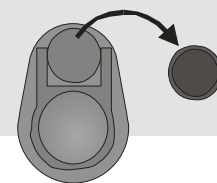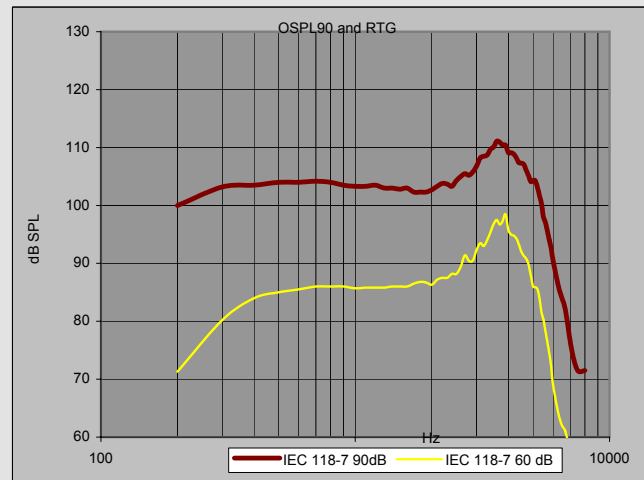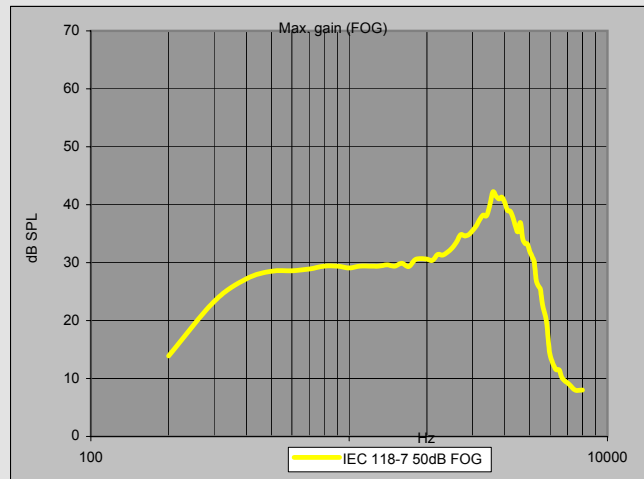

Before programming, please remove the device cover plate

Programming:

Bruckhoff fitting software and Siemens programming cable (CONNEXX cable grey)

Technical change without prior notice